## Bedienungsanleitung Kat. Nr. 50.1010

### BIG STEP - Personenwaage

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
## Instruction manual Kat. Nr. 50.1010
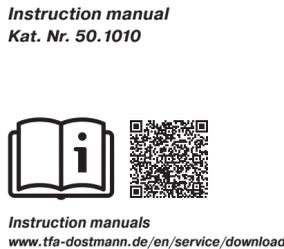
### BIG STEP - Bathroom scale

Thank you for choosing this instrument from the firm TFA.

#### 1. Before you use this product

- **Please make sure you read the instruction manual carefully.** This information will help you to familiarise yourself with your new device, to learn all of its functions and parts, to find out important details about its first use, how to operate it and how to get advice in the event of a breakdown.
- **Following and respecting the instructions in your manual will prevent damage to your instrument and loss of your statutory rights arising from defects due to incorrect use.**
- **We shall not be liable for any damage occurring as a result of non-following of these instructions. Likewise, we take no responsibility for any incorrect readings or for any consequences resulting from them.**
- **Please take particular note of the safety advice!**
- **Please keep this instruction manual safe for future reference.**

#### 2. Delivery contents

- Bathroom scale
- Batteries 3 x 1,5 V AAA
- Instruction manual

#### 3. Field of operation and all the benefits of your new instrument at a glance

- XXL scale: Extra wide step surface and high load capacity
- Automatic power-on function – no tapping required
- Automatic power-off function
- Modern flat design
- Tempered safety-glass surface
- High precision weighing technology (strain gauge sensors)
- Backlight

#### 4. For your safety

- This product is exclusively intended for the field of application described above. It should only be used as described within these instructions.
- Unauthorised repairs, modifications or changes to the product are prohibited.
- Before starting a treatment or diet discuss with your doctor first.

#### **⚠ Caution!** Risk of injury:

- Do not stand on the scale with wet feet and especially not when the surface of the scale is wet. There is a dangerous risk of slipping.
- Do not step on the edge of the scale: Danger of tilting!
- Keep the product and the batteries out of reach of children.
- Batteries contain harmful acids and may be hazardous if swallowed. If a battery is swallowed, this can lead to serious internal burns and death within two hours. If you suspect a battery could have been swallowed or otherwise caught in the body, seek medical help immediately.
- Batteries must not be thrown into a fire, short-circuited, taken apart or recharged. Risk of explosion!
- Low batteries should be changed as soon as possible to prevent damage caused by leaking. Never use a combination of old and new batteries together, nor batteries of different types.
- Wear chemical-resistant protective gloves and safety glasses when handling leaking batteries.


#### **⚠ Important information on product safety!**

- Do not place your product near extreme temperatures, vibrations or shocks.
- Do not place any objects on the unit.
- Protect it from humidity.

#### 5. Set up and operation

- Open the battery compartment at the back and insert the three batteries 1,5 V AAA with the polarity as indicated.
- All LCD segments will light up. "0.0 kg" (default: kg) appears on the display.
- Your scale is now ready to use.
- Press the UNIT button at the back and select your desired weight unit.
- Place the scale on a flat solid surface.

## Instruction manual Kat. Nr. 50.1010
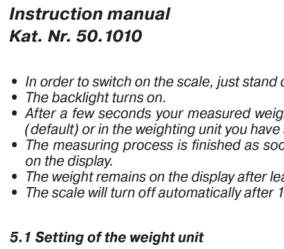
### BIG STEP - Balance pèse-personne

Thank you for choosing this instrument from the firm TFA.

#### 1. Before you use this product

- **Please make sure you read the instruction manual carefully.** This information will help you to familiarise yourself with your new device, to learn all of its functions and parts, to find out important details about its first use, how to operate it and how to get advice in the event of a breakdown.
- **Following and respecting the instructions in your manual will prevent damage to your instrument and loss of your statutory rights arising from defects due to incorrect use.**
- **We shall not be liable for any damage occurring as a result of non-following of these instructions. Likewise, we take no responsibility for any incorrect readings or for any consequences resulting from them.**
- **Please take particular note of the safety advice!**
- **Please keep this instruction manual safe for future reference.**

#### 2. Delivery contents

- Bathroom scale
- Batteries 3 x 1,5 V AAA
- Instruction manual

#### 3. Field of operation and all the benefits of your new instrument at a glance

- XXL scale: Extra wide step surface and high load capacity
- Automatic power-on function – no tapping required
- Automatic power-off function
- Modern flat design
- Tempered safety-glass surface
- High precision weighing technology (strain gauge sensors)
- Backlight

#### 4. For your safety

- This product is exclusively intended for the field of application described above. It should only be used as described within these instructions.
- Unauthorised repairs, modifications or changes to the product are prohibited.
- Before starting a treatment or diet discuss with your doctor first.

#### **⚠ Caution!** Risk of injury:

- Do not stand on the scale with wet feet and especially not when the surface of the scale is wet. There is a dangerous risk of slipping.
- Do not step on the edge of the scale: Danger of tilting!
- Keep the product and the batteries out of reach of children.
- Batteries contain harmful acids and may be hazardous if swallowed. If a battery is swallowed, this can lead to serious internal burns and death within two hours. If you suspect a battery could have been swallowed or otherwise caught in the body, seek medical help immediately.
- Batteries must not be thrown into a fire, short-circuited, taken apart or recharged. Risk of explosion!
- Low batteries should be changed as soon as possible to prevent damage caused by leaking. Never use a combination of old and new batteries together, nor batteries of different types.
- Wear chemical-resistant protective gloves and safety glasses when handling leaking batteries.

#### 9. Specifications

Measuring range:	6 – 200 kg / 440 lb / 31 st
Resolution:	0.1 kg / 0.2 lb
Power consumption:	Batteries 3 x 1.5 V AAA (included)
Dimensions:	380 x 27 x 330 mm
Weight:	2762 g (instrument only)

TFA Dostmann GmbH & Co.KG, Zum Ottersberg 12, D-97877 Wertheim, Germany  
No part of this manual may be reproduced without written consent of TFA Dostmann. The technical data are correct at the time of going to print and may change without prior notice. The latest technical data and information about this product can be found in our homepage by simply entering the product number in the search box.  
[www.tfa-dostmann.de](http://www.tfa-dostmann.de)

03/21

## Mode d'emploi Kat. Nr. 50.1010

### BIG STEP - Balance pèse-personne

Thank you for choosing this instrument from the firm TFA.

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