**TFA**

CE 0197

D

## Digitales Fieberthermometer MT-118

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Lot-Nummer/Chargennummer

Hersteller: Vega Technologies Inc., Yangwu District, Dalang Town, Dongguan City, Guangdong, China

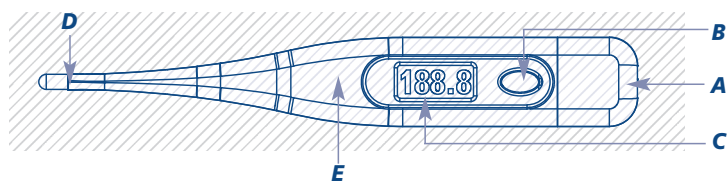
EU-Beauftragter: Biolife Ltd., No.332, Kifisias Av., 15235 Halandri, Athens, Greece

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www.tfa-dostmann.de  
E-Mail: info@tfa-dostmann.de  
Made in P.R.C.

Revision Date: 2021-04-21 D/GB

**TFA®**

0197

## Digital medical thermometer MT-118

Thank you for choosing this instrument from TFA.

### 1. Before you use this product

- Please make sure you read the instruction manual carefully.
- Following and respecting the instructions in your manual will prevent damage to your instrument and loss of your statutory rights arising from defects due to incorrect use.
- We shall not be liable for any damage occurring as a result of non-following of these instructions.
- Likewise, we take no responsibility for any incorrect readings or for any consequences resulting from them.
- Please take particular note of the safety advice!
- Please keep this instruction manual safe for future reference.

### 2. For your safety

- This device is only to be used for measuring human body temperature through oral, rectal or axillary. Do not attempt to take temperatures at other sites, such as in the ear, as it may result in false readings and may lead to injury.
- Before each measurement check intactness of tip and display. Use the unit only when intact.
- Unauthorised repairs, modifications or changes to the product are prohibited.
- Operation in strong electromagnetic fields, such as next to a mobile phone, can lead to malfunction. Keep a minimum distance of 3.3 m from such devices when using this device.
- The minimum measurement time until the beep must be maintained without exception.
- Consider that different measurement locations may require continued measuring even after the beep, see section «Measuring methods / Normal body temperature».
- Do not attempt rectal measurements on persons with rectal disorders. Doing so may aggravate or worsen the disorder.
- **WARNING:** The measurement result given by this device is not a diagnosis! Do not rely on the measurement result only. Discuss the measured values with your doctor.
- The manufacturer recommends to verify the accuracy by an authorised laboratory every 2 years.

**Caution!**  
Risk of injury:

- The instrument contains small parts that could be swallowed by children. Keep this instrument and the battery out of reach of children. Use under parental supervision.
- Batteries contain harmful acids and may be hazardous if swallowed. If a battery is swallowed, this can lead to serious internal burns and death within two hours. If you suspect a battery could have been swallowed or otherwise caught in the body, seek medical help immediately.
- Batteries must not be thrown into a fire, short-circuited, taken apart or recharged. Risk of explosion!
- Low batteries should be changed as soon as possible to prevent damage caused by leaking.
- Wear chemical-resistant protective gloves and safety glasses when handling leaking batteries.
- Remove the battery if you do not use the product for a long period of time.

**Important information on product safety!**

- Protect the device from impact and dropping!
- Do not bend the instrument.
- Avoid ambient temperatures above 60 °C. NEVER boil this device!
- Use only the commercial disinfectants listed in the section «Cleaning and Disinfecting» to clean the device to avoid damage to the device.
- The instrument is not waterproof. Protect from moisture.

### 3. Elements and buttons

- A: Battery compartment  
 B: ON/OFF button  
 C: Display for the temperature  
 D: Measuring sensor / measuring tip  
 E: Cleaning and disinfecting area (entire thermometer surface)

### 4. Use

- Disinfect the tip before each use!
- Press and hold the ON/OFF button for 3 seconds.
- When the unit is switched on, a beep sounds and all segments are briefly displayed. The device performs a self-test.
- The last measured temperature (M) appears briefly on the display, then 36,5 °C.
- If you press the ON/OFF button only briefly when switching on, the last measured temperature (M) will not appear.
- L remains on the display and °C flashes. If the room temperature is higher than 32 °C, the room temperature would be displayed instead of L °C. The device is now ready for use.
- If the temperature remains below 32 °C, °C flashes for another 16 seconds, then the alarm signal sounds.
- The actually measured maximum temperature is permanently shown during the measurement and °C flashes.
- If the maximum temperature remains stable, °C stops flashing and an alarm will sound:
  - a) Normal alarm (beep ... beep ... beep ...), body temperature below 37,8 °C
  - b) Fever alarm (beep pip beep ... beep beep beep ...), body temperature 37,8 °C or above
- For axillary measurements, the recommended measuring time has to be observed despite the beep tone.
- Remove the thermometer carefully. Press the ON/OFF button to switch off the instrument. The instrument will turn off automatically after approx. 9 minutes.

#### 4.1 Measuring methods / Normal body temperature

##### In the mouth (oral) / 35.5 - 37.5 °C

- Do not eat or drink anything hot or cold 10 minutes before the measurement. The mouth should remain closed up to 2 minutes before starting a reading.
- Position the thermometer in one of the two pockets under the tongue, to the left or right of the root of the tongue. The measuring sensor (see point D in the drawing) must be in good contact with the tissue. Close your mouth and breathe evenly through the nose to prevent the measurement from being influenced by inhaled/exhaled air.
- Approx. measuring time: 20 seconds!**

##### In the anus (rectal) / 36.6 - 38.0 °C

- Carefully insert the measuring sensor (see point D in the drawing) 2 to 3 cm into the anal aperture. You may consider using a lubricant. If you are unsure of this measurement method, you should consult a professional for guidance/training.
- Approx. measuring time: 10 seconds!**

##### In the armpit (axillary) / 34.7 - 37.3 °C

- To receive more reliable results, we recommend measuring temperature orally or rectally. A minimum measurement time of 5 minutes is recommended regardless of the beep sound for axillary measurement.

### 5. Cleaning and disinfecting

- For disinfection in home use environment, use a 70% Isopropyl alcohol swab or a cotton tissue moistened with 70% Isopropyl alcohol to wipe surface pollutants off the thermometer probe. Always start wiping from the end of the thermometer probe (approx. at the middle of the thermometer) towards the thermometer tip.
- For professional use: You may apply probe covers (available as accessory). Contact your retailer for further details concerning suitable disinfectants for your thermometer model.

## Instruction manual

Instruction manuals

[www.tfa-dostmann.de/en/service/downloads/instruction-manuals](http://www.tfa-dostmann.de/en/service/downloads/instruction-manuals)
**TFA®**
**Art.-No.15.2008**

## Digital medical thermometer MT-118

### 5.1 Battery replacement

- Change the battery when the display starts to fade or the battery symbol appears on the display.
- Remove the battery cover by pulling it backwards.
- Carefully remove the battery with a pointed object and insert a new LR 41 battery. Make sure the polarity is correct (+ pole above).
- Do not pull out the battery holder. Otherwise, the device may be damaged.
- Slide back the battery cover. Pay attention to the sealing ring.

### 6. Troubleshooting

Problem	Solution
No display	→ Switch on the instrument → Ensure battery polarity is correct (+ pole above) → Change the battery
Incorrect display	→ Change the battery → Check the correct position of the thermometer (see: Measurement methods)
Display E	→ Remove the battery. Wait one minute and then reinsert the battery

If your device fails to work despite these measures, contact the retailer where you purchased the product.

### 7. Legal basis for placing of the product in the EU market

- The Digital Thermometer MT-118, Art.-No. 15.2008 is a medical device class 2a according to the EU Medical Device Directive 93/42/EEG. The Digital Thermometer meets the requirements of the standards EN 12470-3:2000+A1:2009, EN 60601-1:2006+A1:2013, EN 60601-1-11:2015, EN 60601-1-2:2015, EN 60601-1-6:2010+A1:2015, EN ISO 15223-1:2016, EN ISO 10993-1,-5:2009, EN ISO 10993-10:2010, EN ISO 14971:2012.

### 8. Waste disposal

This product has been manufactured using high-grade materials and components which can be recycled and reused.

Never dispose of empty batteries and rechargeable batteries in household waste. As a consumer, you are legally required to take them to your retail store or to appropriate collection sites depending on national or local regulations in order to protect the environment.

The symbols for the heavy metals contained are:  
Cd=cadmium, Pb=lead, Hg=mercury

This instrument is labelled in accordance with the EU Waste Electrical and Electronic Equipment Directive (WEEE). Please do not dispose of this instrument in household waste. The user is obligated to take end-of-life devices to a designated collection point for the disposal of electrical and electronic equipment, in order to ensure environmentally-compatible disposal.

### 9. Specifications

Type:	Maximum thermometer
Measuring range:	32 °C...44 °C
Display L °C	Surrounding temperature lower than 32 °C
Display H °C	Surrounding temperature higher than 44 °C
Accuracy:	±0,1 °C at 35 °C...39 °C otherwise ±0,2 °C
Resolution:	0.1 °C
Storage and transport conditions:	Temperature -10 °C...+60 °C / humidity: 25... 90 % rH
Ambient conditions by use:	Temperature 10 °C...40 °C / humidity: 30...85 % rH
IP protection class:	IP22
Battery:	LR 41 included
Average battery life:	2 years
Housing dimension:	19 x 10 x 124 mm
Weight:	10 g (instrument only)

### 10. Explanation of symbols

0197

Declaration of European conformity, number of the notified body (TÜV Rheinland)

Device complies with protection class BF

Consult instructions for use

Storage and transport conditions:  
Temperature: -10 °C...+60 °C  
Humidity: 25...90 % rH

Ambient conditions by use:  
Temperature: +10 °C...+40 °C  
Humidity: 30...85 % rH

Lot number

Manufacturer: Vega Technologies Inc., Yangwu District, Dalang Town, Dongguan City, Guangdong, China

EU Representative: Biolife Ltd., No.332, Kifisias Av., 15235 Halandri, Athens, Greece

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No part of this manual may be reproduced without written consent of TFA Dostmann. The technical data are correct at the time of going to print and may change without prior notice. The latest technical data and information about this product can be found in our homepage by simply entering the product number in the search box.

www.tfa-dostmann.de  
E-Mail: info@tfa-dostmann.de  
Made in P.R.C.

Revision Date: 2021-04-21 D/GB

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Importer: TFA Dostmann GmbH & Co.KG, Zum Ottersberg 12, D-97877 Wertheim, Germany