Bedienungsanleitung
Instruction manual
Mode d'emploi
Istruzioni per l'uso
Gebruiksaanwijzing
Instrucciones de uso

Instruction manuals
www.tfa-dostmann.de/en/service/downloads/instruction-manuals

Kat. Nr. 50.1019.01

Fig. 1

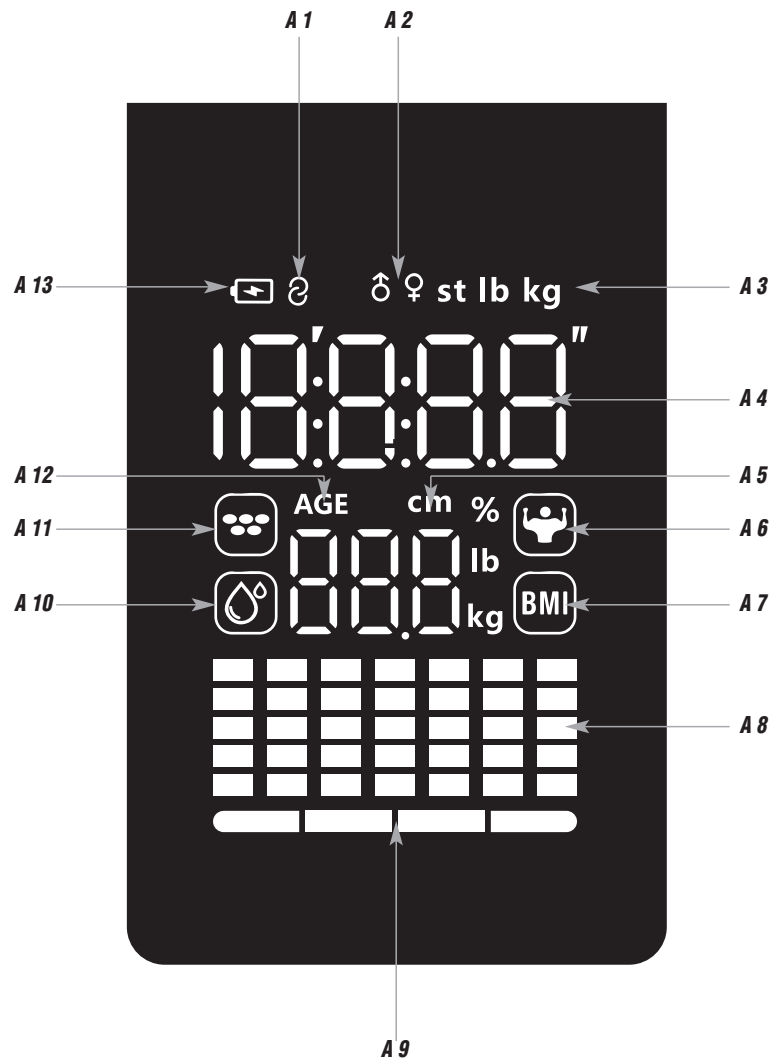
Fig. 2


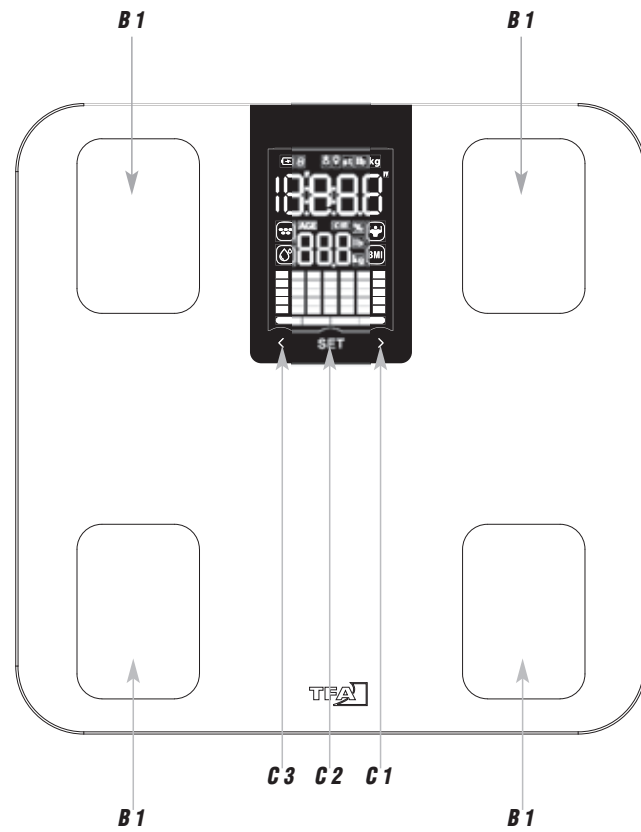
Fig. 3

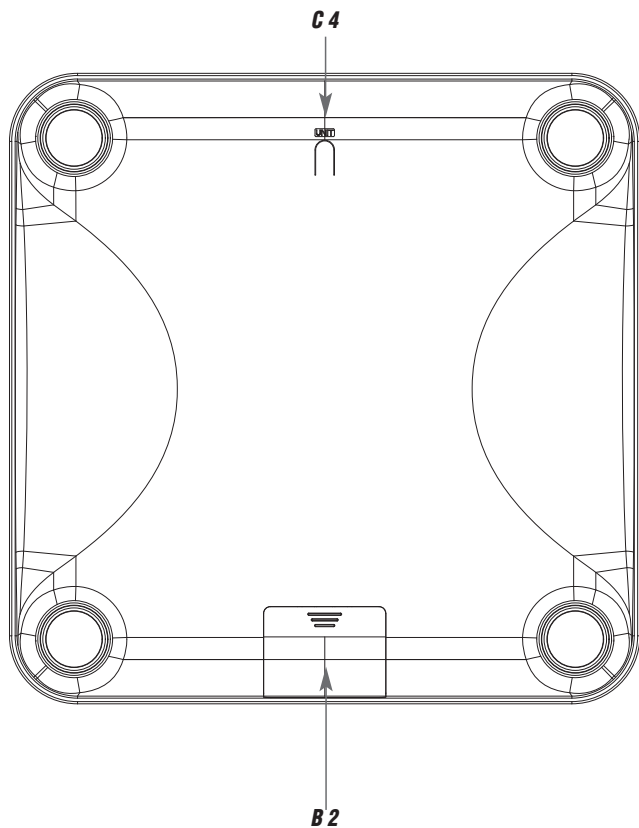Digitale Körper-Analysewaage

D

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Digital body analysis scale

Thank you for choosing this product from TFA Dostmann.

This body composition scale is designed to measure weight and determine body fat, visceral fat, body water, skeletal muscle percentage and muscle mass, bone mass, body age, protein and basal metabolic rate (BMR) using bioelectrical impedance analysis.

A weak alternating current which is harmless to humans is sent through the body by the way of the electrodes placed on the feet and the resistance is then measured. Based on the difference between body components and the various levels of conductivity, the body composition is calculated with the help of the data entered previously, for example, height, sex and age.

A healthy body weight not only depends on the number of kilograms, but also on the muscle mass to fat mass ratio. Having too high a body fat content can increase the risk of cardiovascular diseases. The body analysis scale allows you to regularly check your personal values and promote your health and fitness.

1. Before you use this product

- Please make sure you read the instruction manual carefully .
- The operating instructions are enclosed with the device or can be downloaded at www.tfa-dostmann.de/en/service/downloads/instruction-manuals
- This product should only be used as described within these instructions.
- Following and respecting the instructions in your manual will prevent damage to your instrument and loss of your statutory rights arising from defects due to incorrect use. We shall not be liable for any damage occurring as a result of non-following of these instructions. Likewise, we take no responsibility for any incorrect readings or for any consequences resulting from them.
- Please pay particular attention to the safety notices!
- Please keep this instruction manual safe for future reference.

2. Delivery contents

- Body analysis scale
- Instruction manual

3. Range of application and all the benefits of your new instrument at a glance

- Body composition analysis with many health metrics such as weight, BMI, body fat, visceral fat, body water, skeletal muscle percentage, muscle mass, bone mass, body age, protein and basal metabolic rate (BMR) using bioelectrical impedance analysis
- Data transmission and analysis via app
- Memory storage capacity for four people
- Additional guest user
- Automatic person recognition
- Metabolic rate (BMR) using bioelectrical impedance analysis
- Coloured LED display
- Automatic power on without touching
- Auto switch-off function
- Calibration function
- Tempered safety glass platform

4. For your safety

- All data on body composition are provided for reference purposes only.

Digital body analysis scale

- Before starting a treatment or diet discuss with your doctor first.
- Not suitable for people with electronic implants (pacemakers, etc.).
- Children under 16 and pregnant women should only use the scale to measure their body weight (settings in the app) and should wear shoes during the weighing process.

WARNING!
Risk of injury!

- Do not stand on the scale with wet feet, nor when the surface of the scale is wet. **Risk of slipping!**
- Do not step on the edge of the scale: **Danger of tilting!**
- Keep this device and the batteries out of reach of children.
- Small parts can be swallowed by children (under three years old).
- Batteries contain harmful acids and may be hazardous if swallowed. If you suspect a battery could have been swallowed or otherwise caught in the body, seek medical help immediately.
- Batteries must not be thrown into a fire, short-circuited, taken apart or recharged. **Risk of explosion!**
- Low batteries should be changed as soon as possible to prevent damage caused by leaking. Make sure the polarities are correct. Never use a combination of old and new batteries together, nor batteries of different types. Remove the batteries if the device will not be used for an extended period of time.
- Avoid contact with skin, eyes and mucous membranes when handling leaking batteries. In case of contact, immediately rinse the affected areas with water and consult a doctor.

CAUTION!
Important information on product safety!

- Do not expose the device to extreme temperatures, vibrations or shocks.
- Do not place any objects on the unit.
- Clean the device with a soft damp cloth. Do not use solvents or scouring agents.
- Protect it from moisture!

5. Elements

A: Display (Fig. 1):

- A 1:** Wireless connection symbol
- A 3:** Units
- A 5:** Height
- A 7:** BMI
- A 9:** Body type
- A 11:** Body fat
- A 13:** Battery symbol

- A 2:** Sex
- A 4:** Weight
- A 6:** Muscle mass
- A 8:** Analysis bar
- A 10:** Body water
- A 12:** Age

B: Housing (Fig. 2+3):

- B 1:** Electrodes
- B 2:** Battery compartment

C: Buttons (Fig. 2+3):

- C 1:** > button
- C 2:** SET button
- C 3:** < button
- C 4:** UNIT button

Digital body analysis scale

6. Operation

6.1 Getting started

- Open the battery compartment and insert three new AAA 1.5 V batteries.
- Close the battery compartment.
- Press the **UNIT** button on the back of the scale to switch between kg (default), lb and st as the unit of weight.
- **Note:** If the US or UK weight system is selected, height is displayed in feet (ft).

6.2 Calibration

- The scale must be calibrated when using it for the first time, after changing the batteries, or if it has been moved to a different location.
- Step onto the scale briefly to switch it on, then step off the weighing platform.
- The calibration is finished, as soon as 0.00 kg appears on the display.

6.3 Installing the app and pairing the scale with your smartphone

- Search for "Fitdays" in the App Store or Google Play, or scan the following QR code.

- Install the "Fitdays" app on your smartphone.
- Sign up in the app.
- Enter your personal details and confirm them.
- Activate wireless connection on your smartphone and return to the app.
- Tap the + icon in the top right-hand corner of the app and select 'My Device'.
- Next, tap the '+ Add device' button and follow the instructions in the app.
- **NOTES:**
 - By registering and entering your personal details, you are automatically set up as the key user (U-1).
 - Perform the first measurement with the app connected to pair the device. All subsequent measurements can also be carried out offline without connection to the app.

7. Using the scale with the analysis function

7.1 To ensure reliable measurements

- Always stand on the scale barefoot. Feet must be firmly positioned on the electrodes.
- If possible, always weigh yourself at the same time of day, after going to the toilet and when not wearing any clothes.
- The measurement very much depends on the user's water balance. Saunas, physical activity, drinking, eating, illnesses and the female menstrual cycle can all impact on the results. Wait a few hours after strenuous exercise.
- **To achieve comparable results, you should always weigh yourself when at the same level of hydration if possible.**

Digital body analysis scale

7.2 Performing the body analysis

- Place the scale on a flat solid surface.
- Stand barefoot on the scale. The display turns on and the scale recognises the user.
- **Note:** If the values of several people with the same measurements are stored, you are prompted to select the right user with the < or > button at the scale or select the user in the app.
- After a few seconds your weight will appear.
- Stand still until the analysis is complete.
- The display will now show your height and the analysis bars will appear. After a few seconds appear:

Body fat percentage

Muscle mass percentage

Body water percentage

BMI (Body-Mass-Index)

- The values are displayed one after the other for three seconds.
- The measured values will be transmitted.
- Open the app to view your detailed readings, track your progress and access additional information.

7.2.1 Interpreting the measured values

- **Please note:** the measured values calculated by the body analysis scale are not precise medical values. Contact a specialist for an exact calculation.
- Fat, muscle and bone mass also consist partly of water, and the percentages may therefore add up to more than 100%.

7.2.2 Analysis bar

- Up to seven analysis bars can be displayed on the screen.
- One analysis bar corresponds to one measurement.
- The height of an analysis bar indicates the change from the previous readings.

7.3 Measuring procedure without body analysis

- Once pregnancy mode is activated in the app, or if you step onto the scales whilst wearing shoes, only your weight will be measured and no body composition analysis will be carried out.
- Place the scale on a flat solid surface.
- To switch on the scale, simply step on the platform.
- After a few seconds your measured weight will be displayed in kg (default) or in the weighing unit you have set.
- The weight measuring is finished as soon as your weight will be flashing in the display.
- The weight remains displayed after you step off the scale.
- The scale will turn off automatically after approximately 10 seconds.

7.4 Add user

7.4.1 In the app

- Open the app and go to your account.

Digital body analysis scale

- Select the 'Add User' button.
- Please enter the following details:
 - Nickname
 - Birthday
 - Height
 - Target weight
 - Mode
- The registered user information is transmitted to the scale.

7.4.2 At the scale

- Press the **SET** button. The scale will turn on automatically.
- U-1 flashes and the age and sex are shown on the display.
- Choose with the **<** or **>** button the user (U1...U4).
- Press and hold the **SET** button for three seconds to enter the setting mode for the selected user.
- The age 18 (default) flashes on the display. Set your age with the **<** or **>** button.
- Confirm the setting with the **SET** button.
- The symbol for a man (default) appears on the display.
- Set your sex with the **<** or **>** button.
- Confirm the setting with the **SET** button.
- 170 cm (default) flashes on the display.
- Set your height with the **<** or **>** button.
- Confirm the setting with the **SET** button.
- You can now take the measurement.

8. Battery replacement

- Change the batteries when the battery symbol or "Lo" appears on the display.
- Open the battery compartment and insert three new AAA 1.5 V batteries, polarity as illustrated.

9. Troubleshooting

Problem	Solution
No display	→ Step on the scale → Ensure the batteries' polarities are correct → Change the batteries
Display "C"	→ Place the scale on a flat solid surface → Step on the scale again
Display "Err"	→ The weight is too high
No body analysis performed	→ Stand barefoot on the scale → Stand still until the end of the analysis → Before you start weighing yourself enter your personal data
Incorrect indication	→ Change the batteries

If your device fails to work despite these measures, please contact the retailer where you purchased the product from for advice.

Digital body analysis scale

10. Waste disposal

Dispose of the packaging in an environmentally friendly manner using the collection systems that have been set up.

Disposal of the electrical device

Remove non-permanently installed batteries and rechargeable batteries from the device and dispose of them separately.

This product is labelled in accordance with the EU Waste Electrical and Electronic Equipment Directive (WEEE).

This product must not be disposed of in ordinary household waste. As a consumer, you are required to take end-of-life devices to a designated collection point for the disposal of electrical and electronic equipment, in order to ensure environmentally-compatible disposal. The return service is free of charge. Observe the current regulations in place!

Please note that you are responsible for deleting any personal data from the old devices you are disposing of.

Disposal of the batteries

Never dispose of empty batteries and rechargeable batteries with ordinary household waste. They contain pollutants which, if improperly disposed of, can harm the environment and human health. As a consumer, you are required by law to take them to your retail store or to an appropriate collection site depending on national or local regulations in order to protect the environment. The return service is free of charge.

The symbols for the contained heavy metals are:

Cd=cadmium, Hg=mercury, Pb=lead.

11. Specifications

Measuring range	5 – 180 kg / 396 lb / 28 st
Resolution	0.05 kg / 0.1 lb
Accuracy:	±300 g @ 5 kg...60 kg ±400 g @ 61 kg...120 kg ±500 g @ 121 kg...180 kg
Operating temperature and humidity	Temperature: 10°C...40°C Humidity: 40...80 %rH
Auto switch-off function	After 10 seconds
Power consumption	3 x AAA 1.5 V batteries (not included)
Housing dimension	278 x 278 x 24 mm
Weight	1274 g (device only)

Digital body analysis scale**(GB)**

No part of this manual may be reproduced without written consent of TFA Dostmann. The technical data are correct at the time of going to print and may change without prior notice.

The latest technical data and information about this product can be found in our homepage by simply entering the product number in the search box.

EU Declaration of conformity

Hereby, TFA Dostmann declares that the radio equipment type 50.1019 is in compliance with Directive 2014/53/EU. The full text of the EU declaration of conformity is available at the following Internet address: www.tfa-dostmann.de/service/downloads/ce

Privacy Notice

The device can optionally be used with the third-party app "Fitdays". The use of the app is voluntary. The processing of personal data within the app is carried out exclusively by the respective app provider. This provider's privacy policy and terms of use apply.

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TFA Dostmann GmbH & Co.KG, Zum Ottersberg 12, 97877 Wertheim, Alemania